


Place 7 way Jade Massage Heads
on the spots for **Sciatic Neuralgia**.

 spot for 7 way Jade Massage Heads

- 1: Sciatic Neuralgia, Waist Pain
- 2: Sciatic Neuralgia, Waist & Thigh Pain
- 3: Sciatic Neuralgia, Waist Pain
- 4: Waist & Calf pain
- 5: Waist Pain, Leg Arthritis

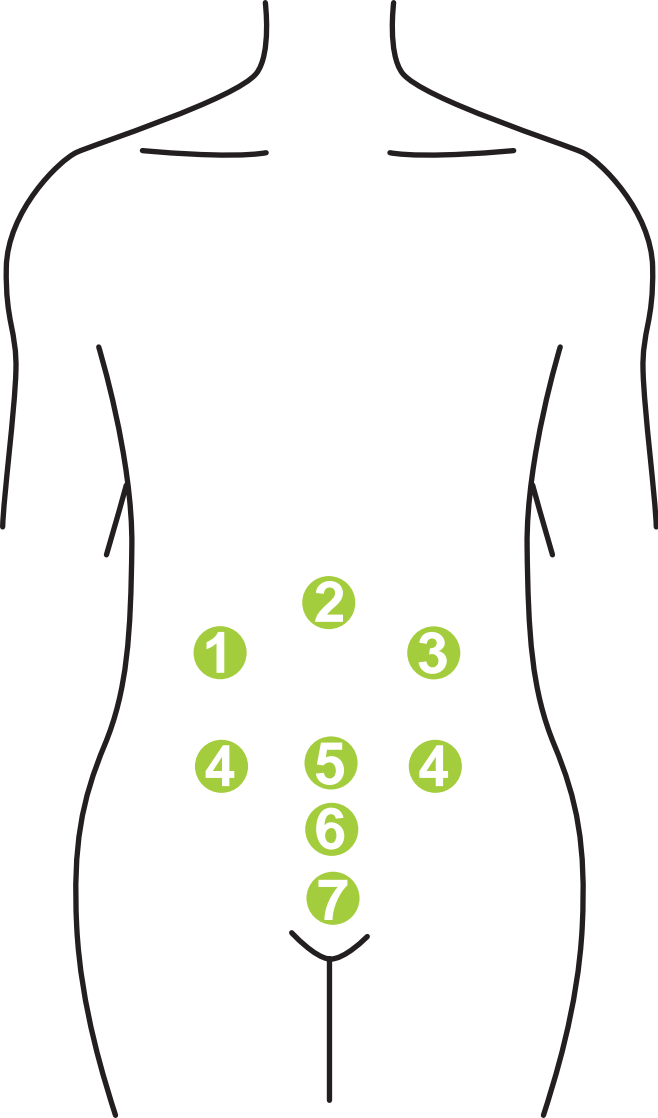

Place 7 way Jade Massage Heads
on the spots for **General Organs.**

● spot for 7 way Jade Massage Heads

1: Liver, Gallbladder

2: Stomach, Spleen


3: Diabetes

4: Constipation

5: Small Intestines

6: Urinary Bladder, Prostate

7: Uterus

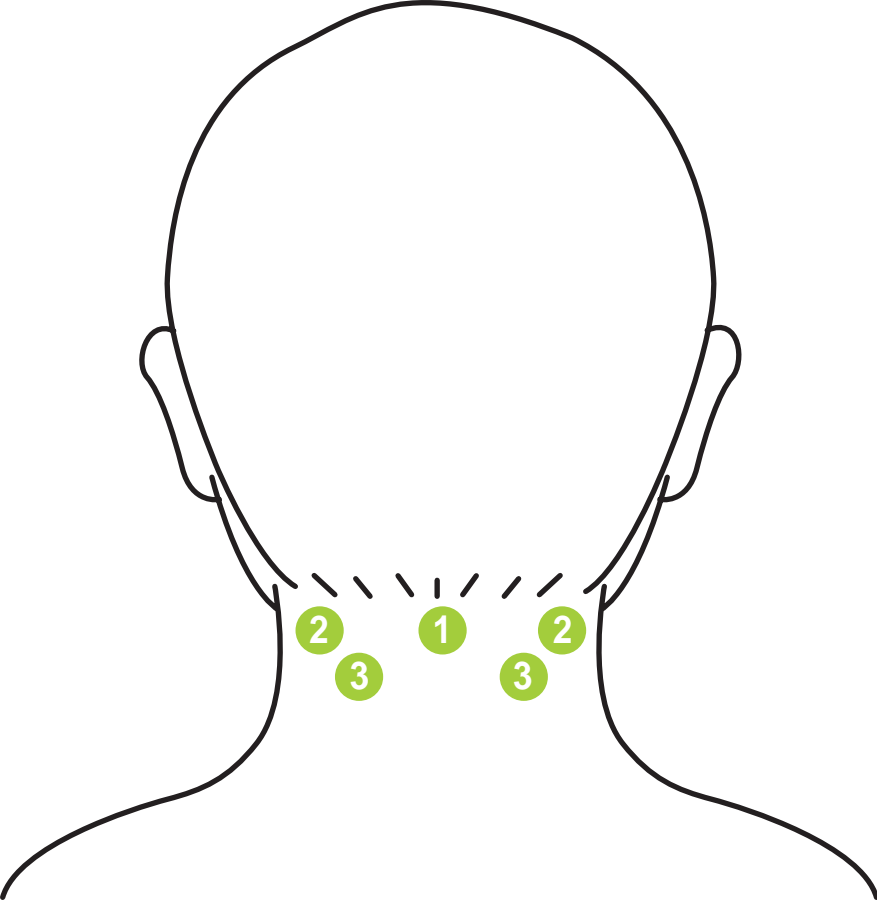


Place 7 way Jade Massage Heads on the spots for **Headache**.



spot for 7 way Jade Massage Heads




- 1: General Headache, Increase Blood Circulation to the Brain**
- 2: Paralysis, Dizziness**
- 3: Headache, Neck Pain**

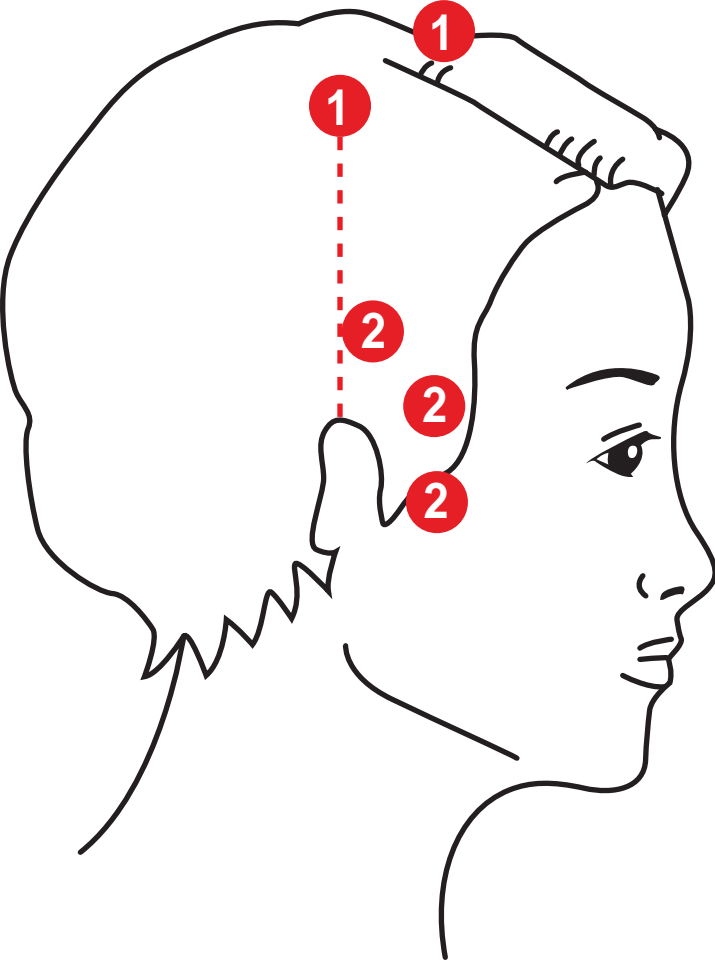

Place 2 way Jade Massage Heads on the spots for **Migraine**.

● spot for 2 way Jade Massage Heads

- 1: Headache,
Cerebral Hemorrhage
(For Position 1, do not set the temperature more than 45 Celsius)
- 2: Migraine

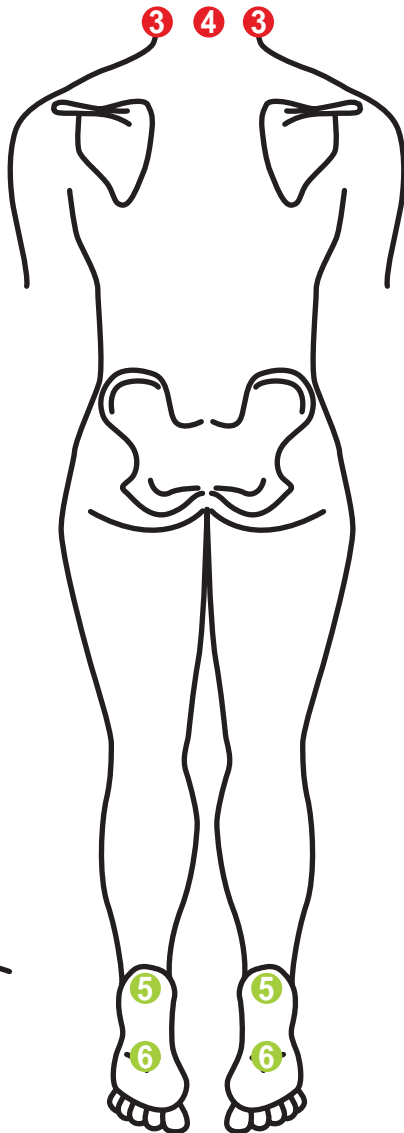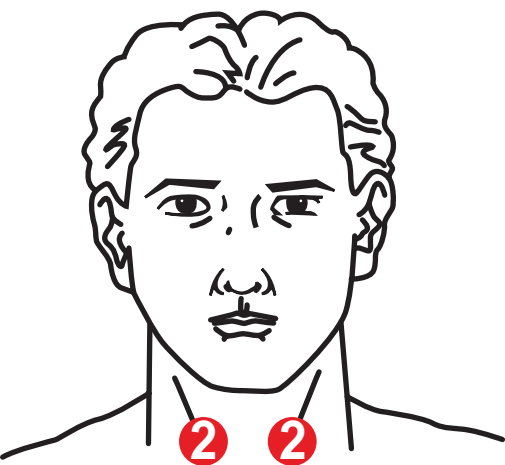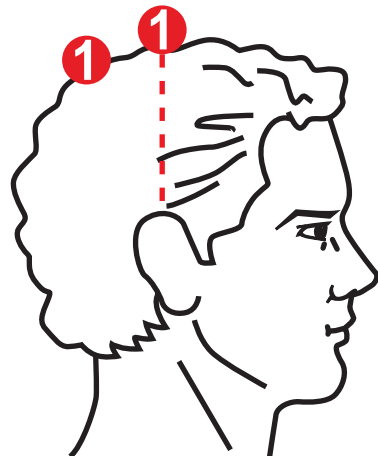


Place 2 & 7 way Jade Massage Heads on the spots for **High Blood Pressure**.

● spot for 2 way Jade Massage Heads

● spot for 7 way Jade Massage Heads

- 1: **Cerebral Hemorrhage, Stroke** (For position 1, Do not set the temperature more than 45 Celsius)
- 2: **High Blood Pressure**
- 3: **Stroke**
- 4: **Increase Blood Circulation to the Brain**
- 5: **Recovering Fatigue**
- 6: **Heart Disease, High Blood Pressure**

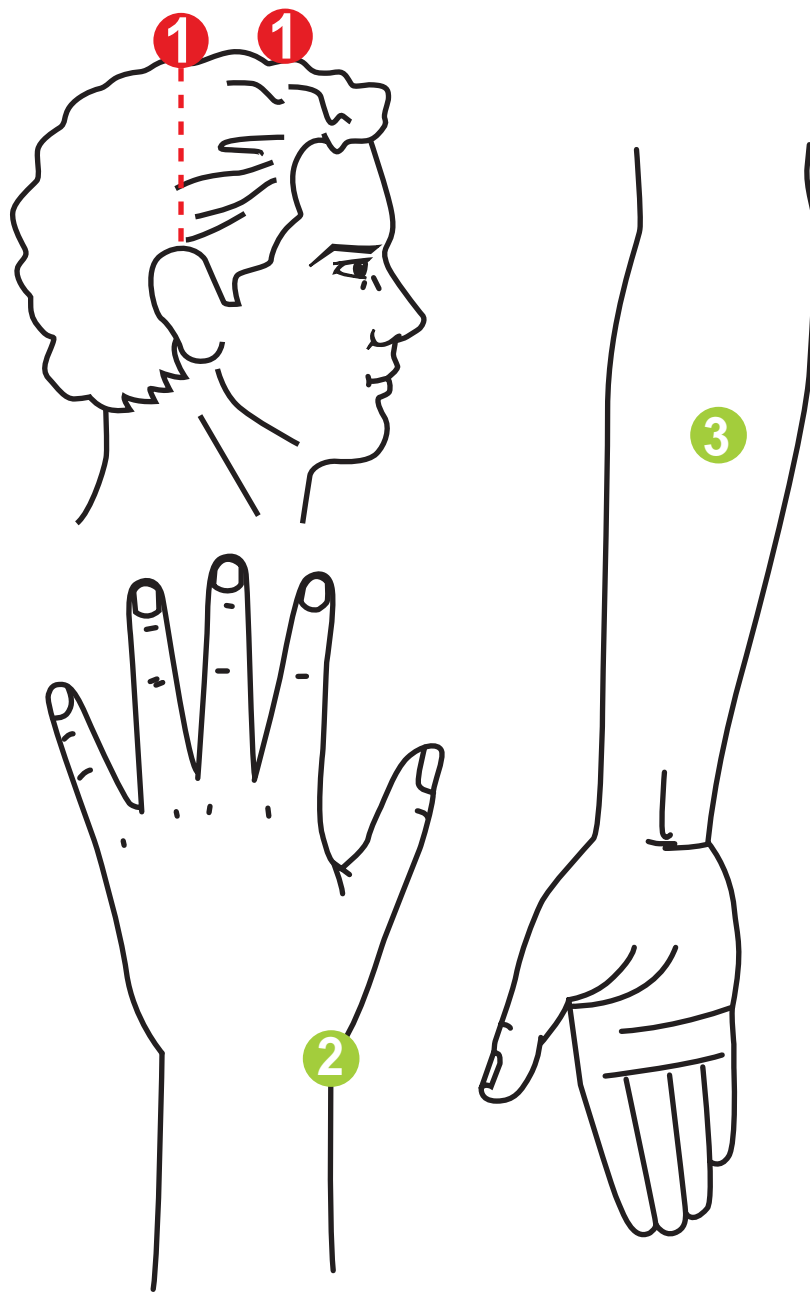

Place 2 & 7 way Jade Massage Heads on the spots for **Low Blood Pressure**.


 spot for 2 way Jade Massage Heads

 spot for 7 way Jade Massage Heads


- 1: Anemia**
(For position 1, Do not set the temperature more than 45 Celsius)
- 2: Cerebral Anemia**
- 3: Improve Blood Circulation**

Place 2 way Jade Massage Heads on the spots for **Eye Problems**.

● spot for 2 way Jade Massage Heads

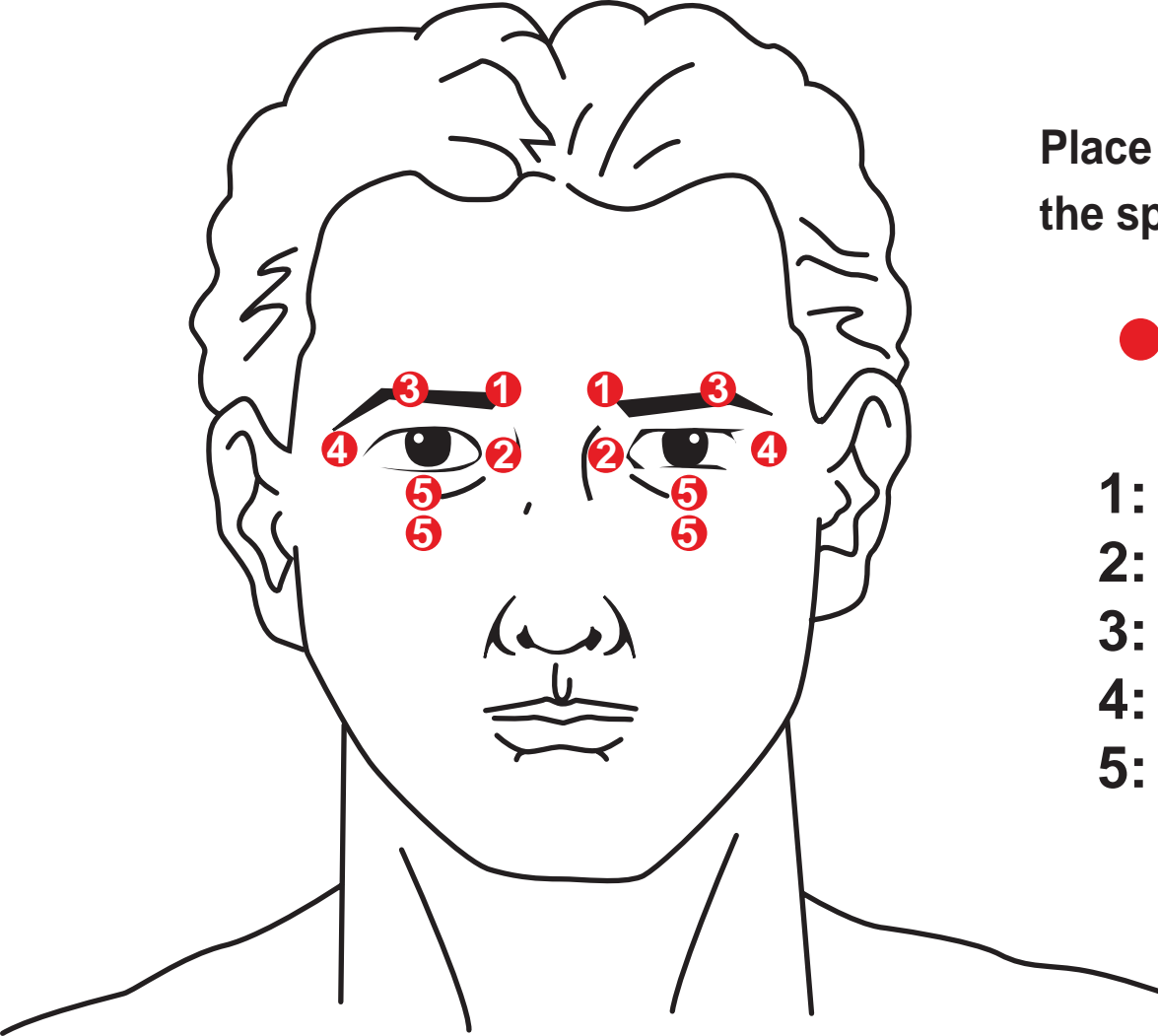


- 1: Tearing Eyes
- 2: Stimulate Lacrimal Gland
- 3: Glaucoma, Sty
- 4: Improve Eyesight
- 5: Conjunctivitis, Improve Eyesight

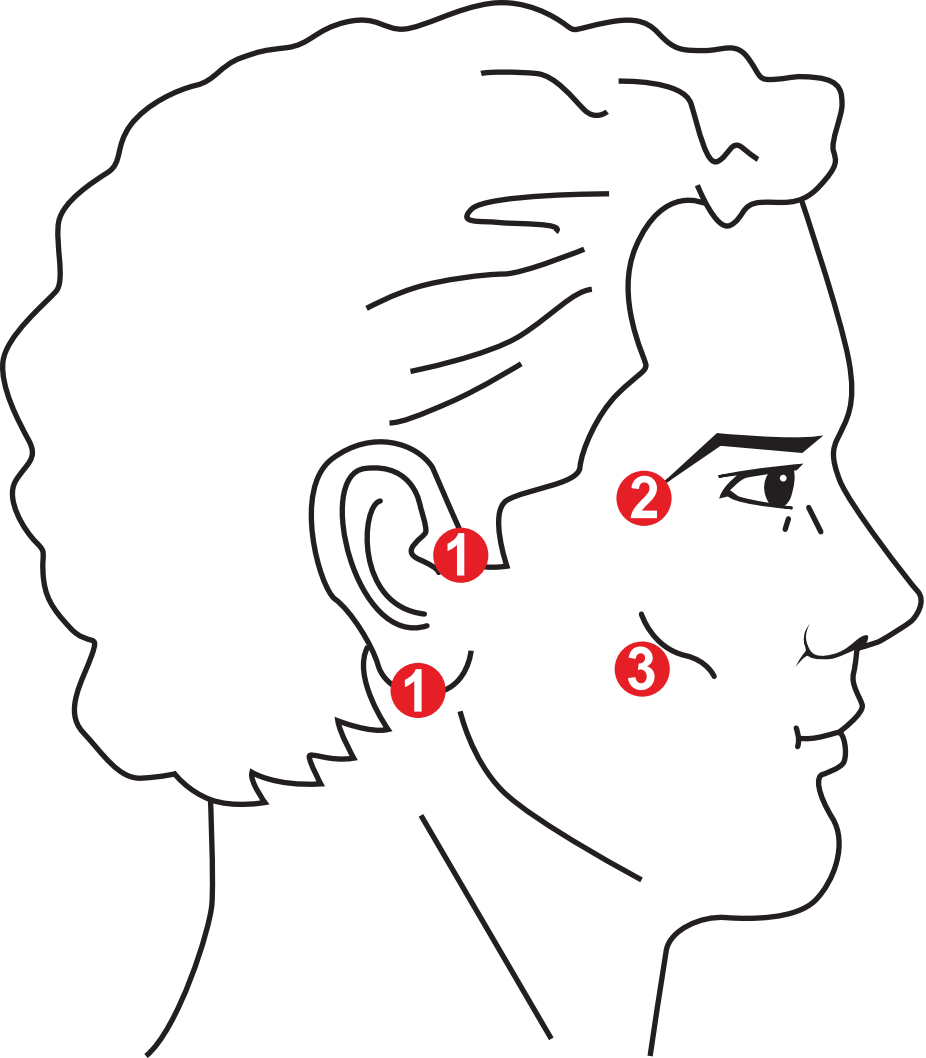

Place 2 way Jade Massage Heads on the spots for **Ear Problems**.

● spot for 2 way Jade Massage Heads


- 1: Tinnitus
Improve Hearing
- 2: Improve Eyesight
- 3: Facial Paralysis

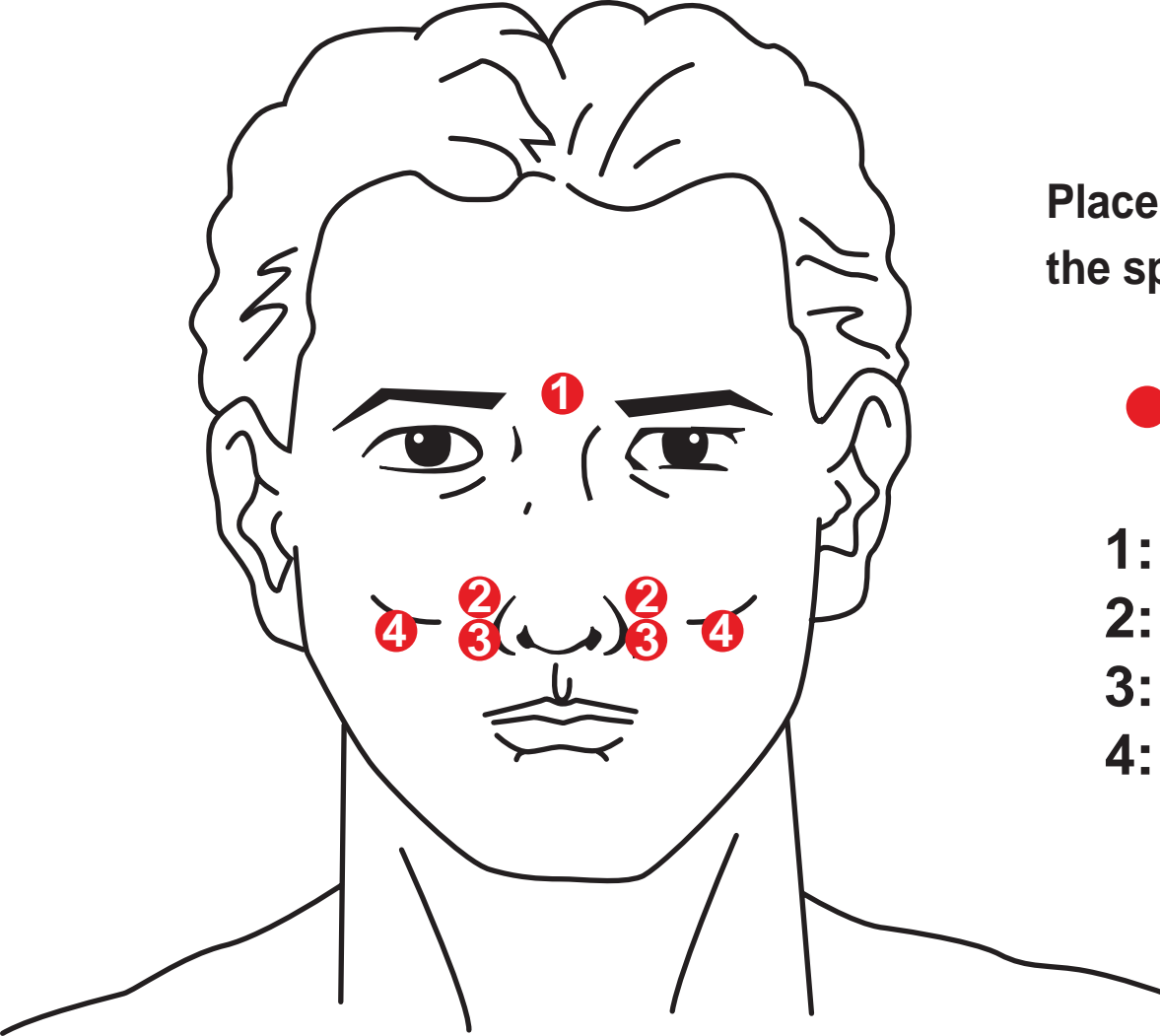

Place 2 way Jade Massage Heads on the spots for **Empyema**.

● spot for 2 way Jade Massage Heads


- 1: Headache
- 2: Allergy
- 3: Empyema
- 4: Facial Paralysis

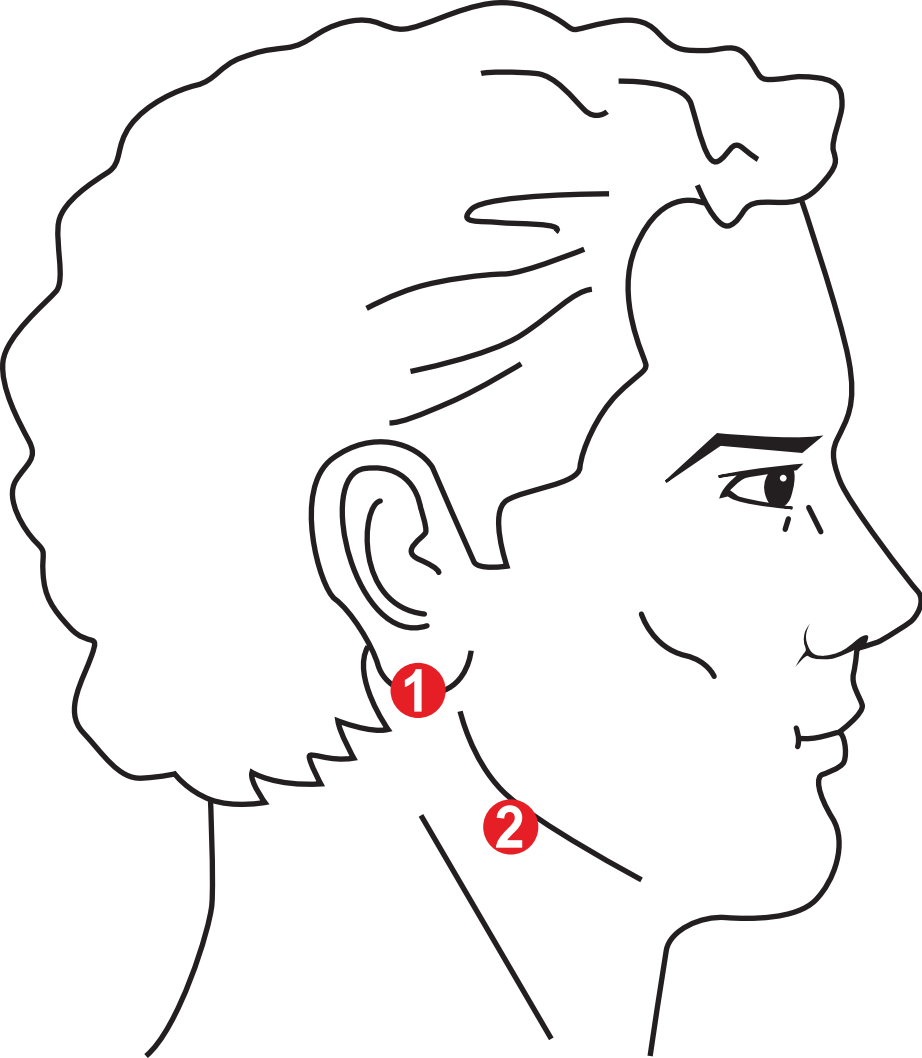

Place 2 way Jade Massage Heads on the spots for **Teeth Problems**.

● spot for 2 way Jade Massage Heads


- 1: Improve Hearing
- 2: Toothache

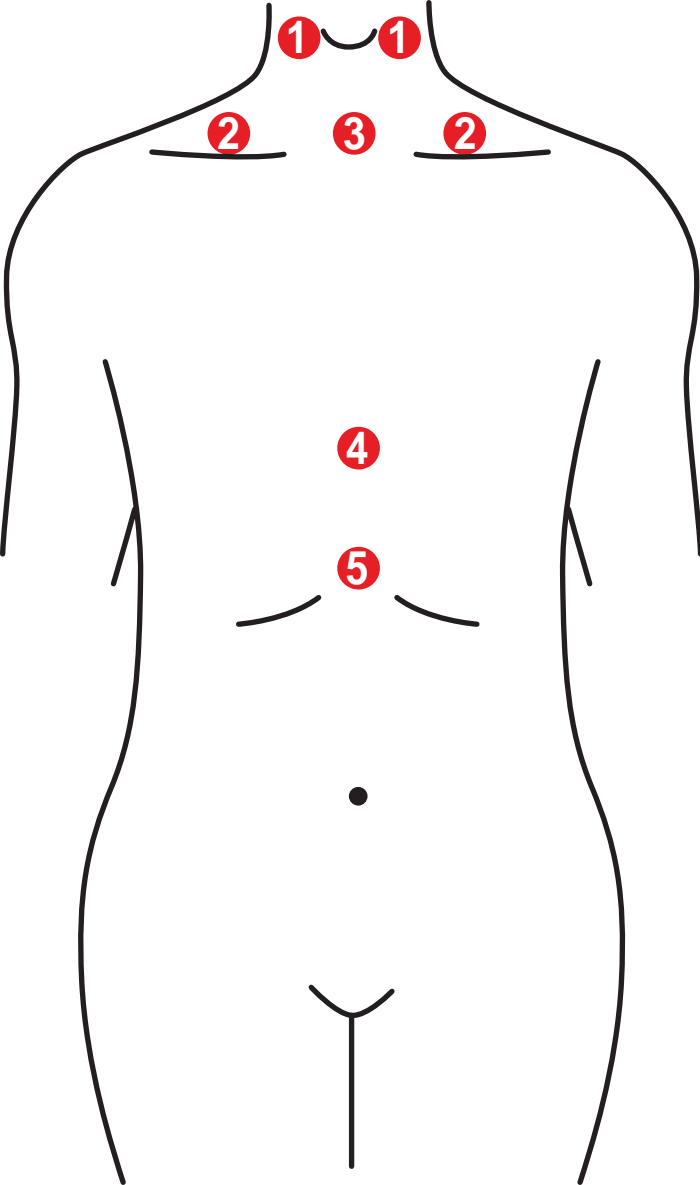

Place 2 way Jade Massage Heads on the spots for **Bronchial Asthma**.

 spot for 2 way Jade Massage Heads


- 1: Asthma, Uvulitis
- 2: Bronchial Problem, Neck Problems
- 3: Cough
- 4: Chest Pain, Tuberculosis
- 5: Cough, Chest Pain

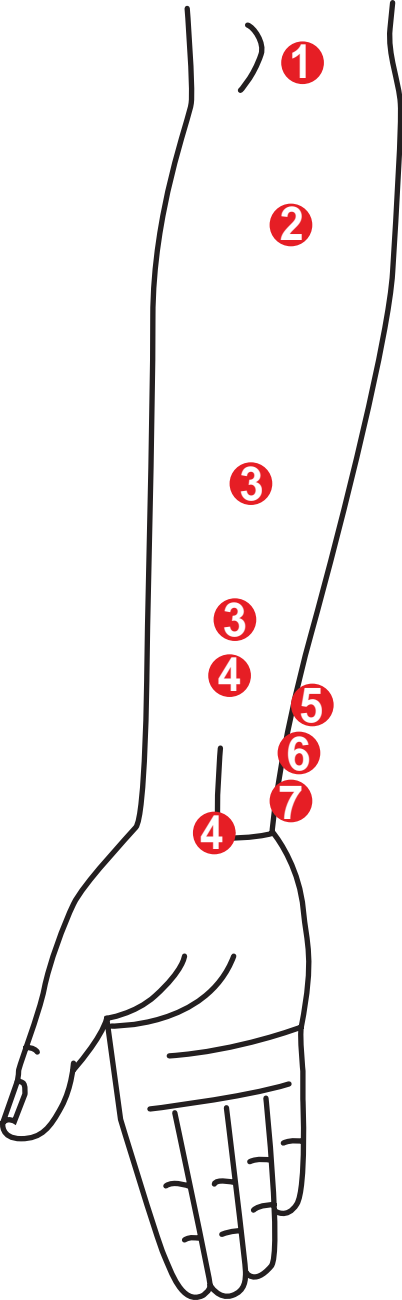


Place 2 way Jade Massage Heads on the spots and 7 way Jade Massage Head on your left arm for **Heart Problems**.

 spot for 2 way Jade Massage Heads

- 1: Pain in the Chest
- 2: Elbow and Finger Problems, Breathing Problems
- 3: Heart Disease, Blood Circulation
- 4: Chest Pain
- 5: Inflammation of Wrist Joints, Chest Pain
- 6: Heart Disease
- 7: Chest Pain

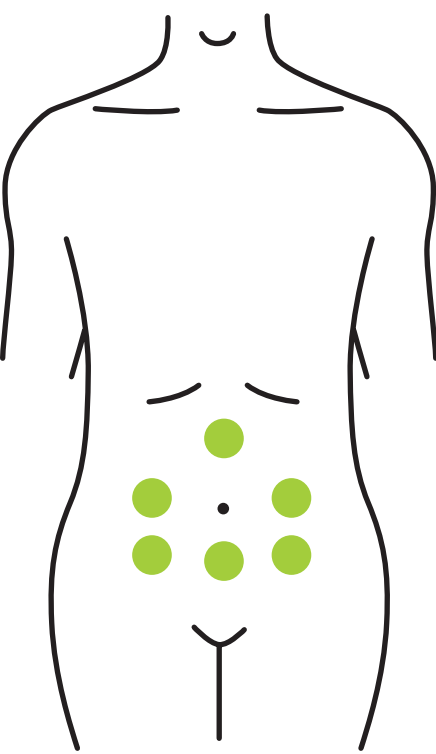

Place 7 way Jade Massage Heads
on the floor and lie on top of it.

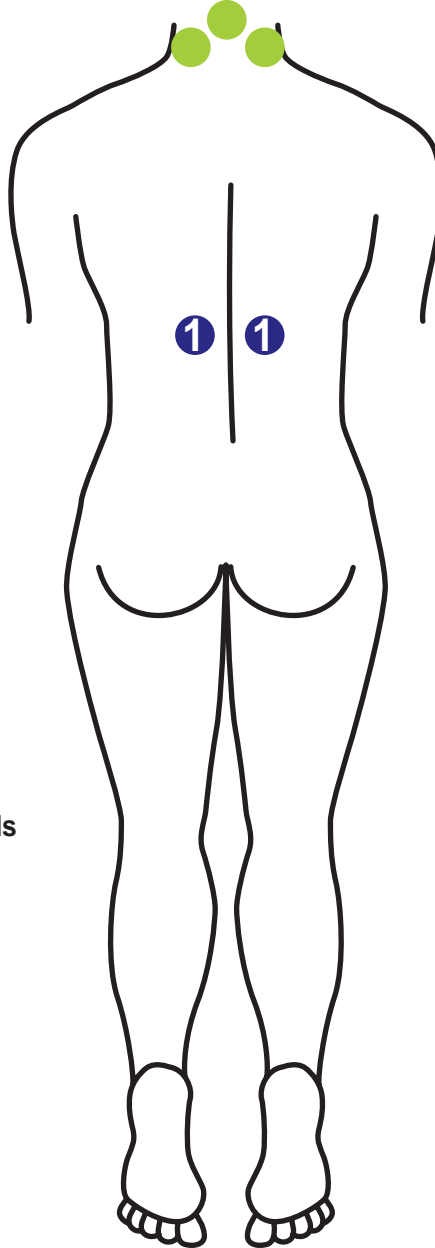


Place 7 & 5 way Jade Massage Heads
on the spots for **Diabetes**.

 spot for 7 way Jade Massage Heads

 spot for 5 way Jade Massage Heads

1: Place 5 way Jade Massage Heads on the
Back, approximately 11th ~12th thoracic
vertebrae.

By placing 5 way Jade Massage Heads on
the above area, You can stimulate the
Pancreas to secrete insulin.

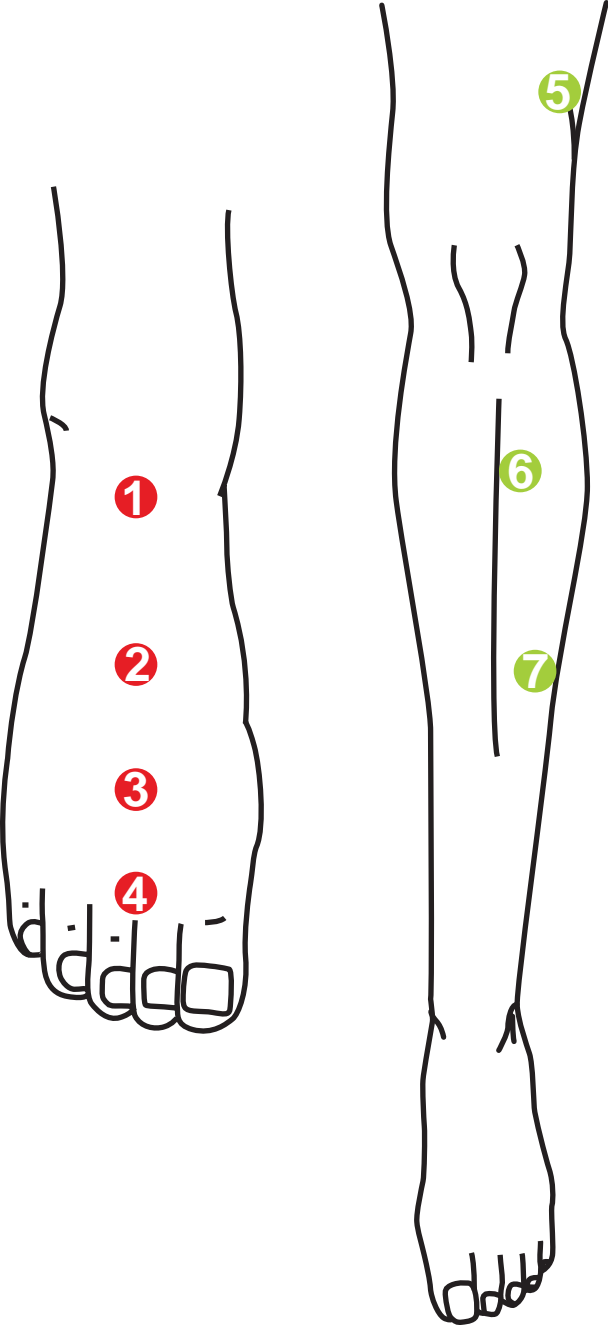


Place 2 & 7 way Jade Massage Heads
on the spots for **Stomach Trouble**.

 spot for 2 way Jade Massage Heads

 spot for 7 way Jade Massage Heads


- 1: Pain in the Ankle, Stomach Trouble
- 2: Intestinal Trouble
- 3: Tonify Stomach & Intestines
- 4: Stomach Trouble
- 5: Knee Arthritis, Release Stomach Stress
- 6: Stomach Trouble
- 7: Release Stomach Stress

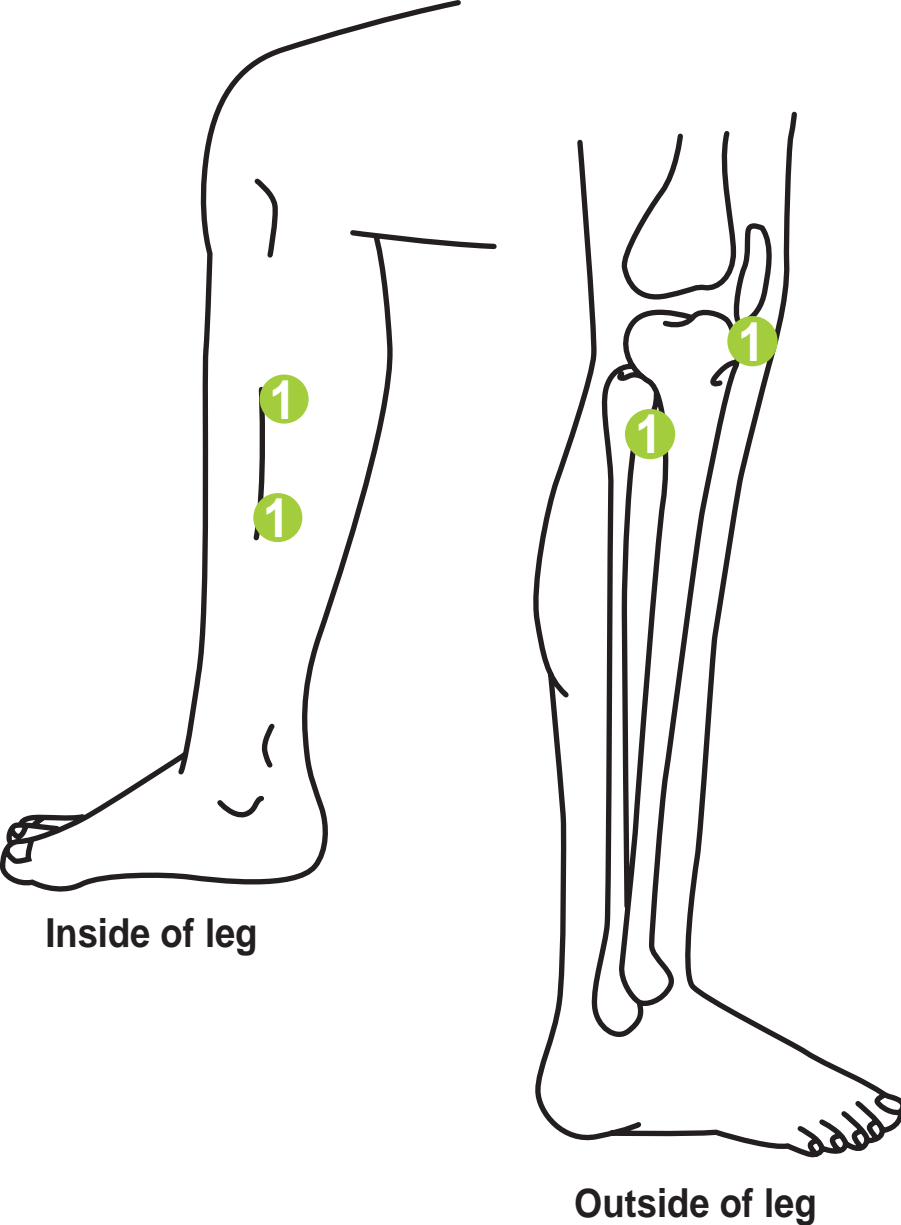

Place 7 way Jade Massage Heads
on the spots for **Liver**.

● spot for 7 way Jade Massage Heads

1: Improve Liver Function

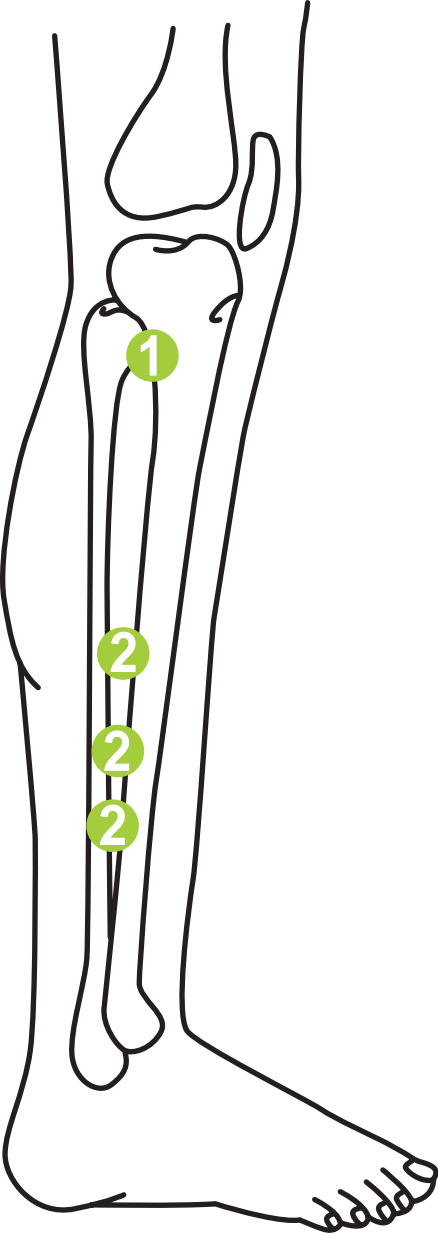


Outside of leg



Place 7 way Jade Massage Heads
on the spots for **Gallbladder**.

● spot for 7 way Jade Massage Heads

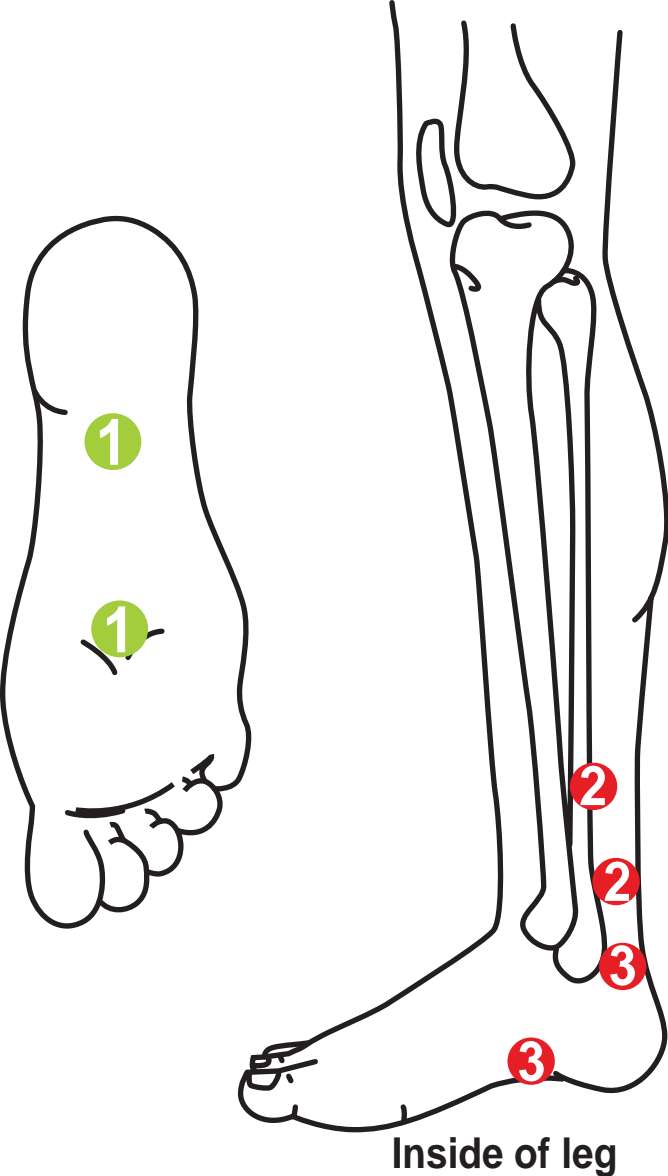
- 1: Improve Liver Function
- 2: Improve Gallbladder Function



Place 2 & 7 way Jade Massage Heads
on the spots for **Kidney Trouble**.

-  spot for 2 way Jade Massage Heads
-  spot for 7 way Jade Massage Heads

- 1: Increase Energy Level of Body**
- 2: Improve Kidney Function**
- 3: Paralysis of the Urinary Bladder**

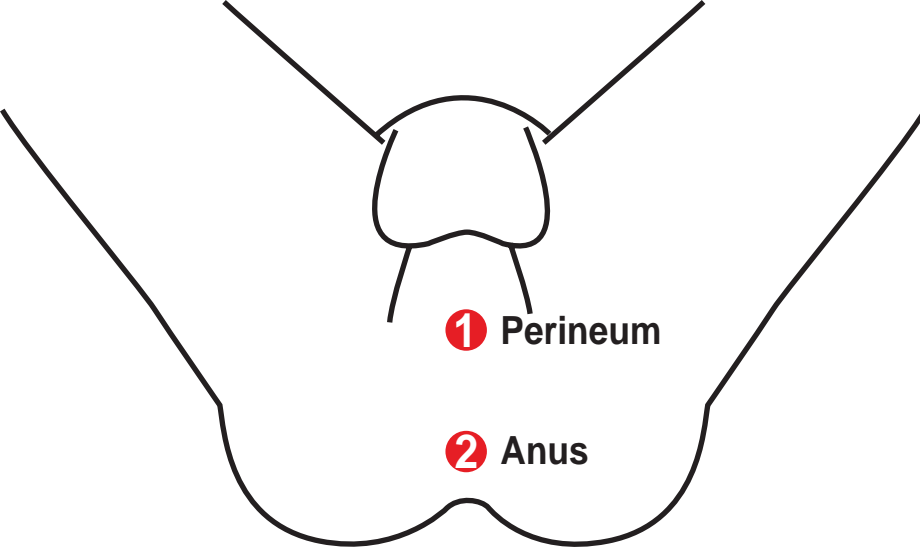


① Perineum

② Anus


Place 2 & 7 way Jade Massage Heads
on the spots for **Prostate Gland**.

● spot for 2 way Jade Massage Heads

● spot for 7 way Jade Massage Heads

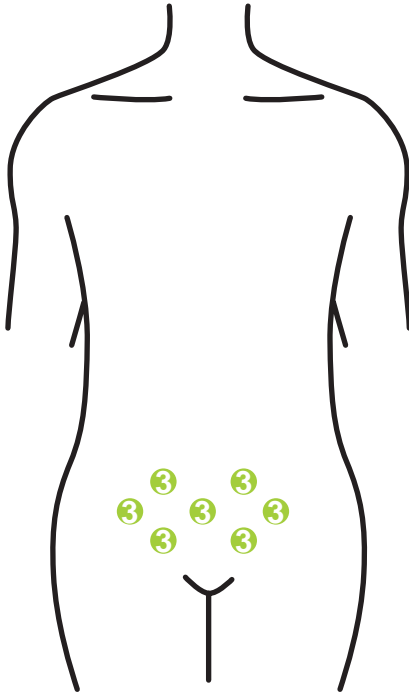


- 1: Prostatomegaly
- 2: Hemorrhoids (Piles)
- 3: Improve Urinary Bladder Function

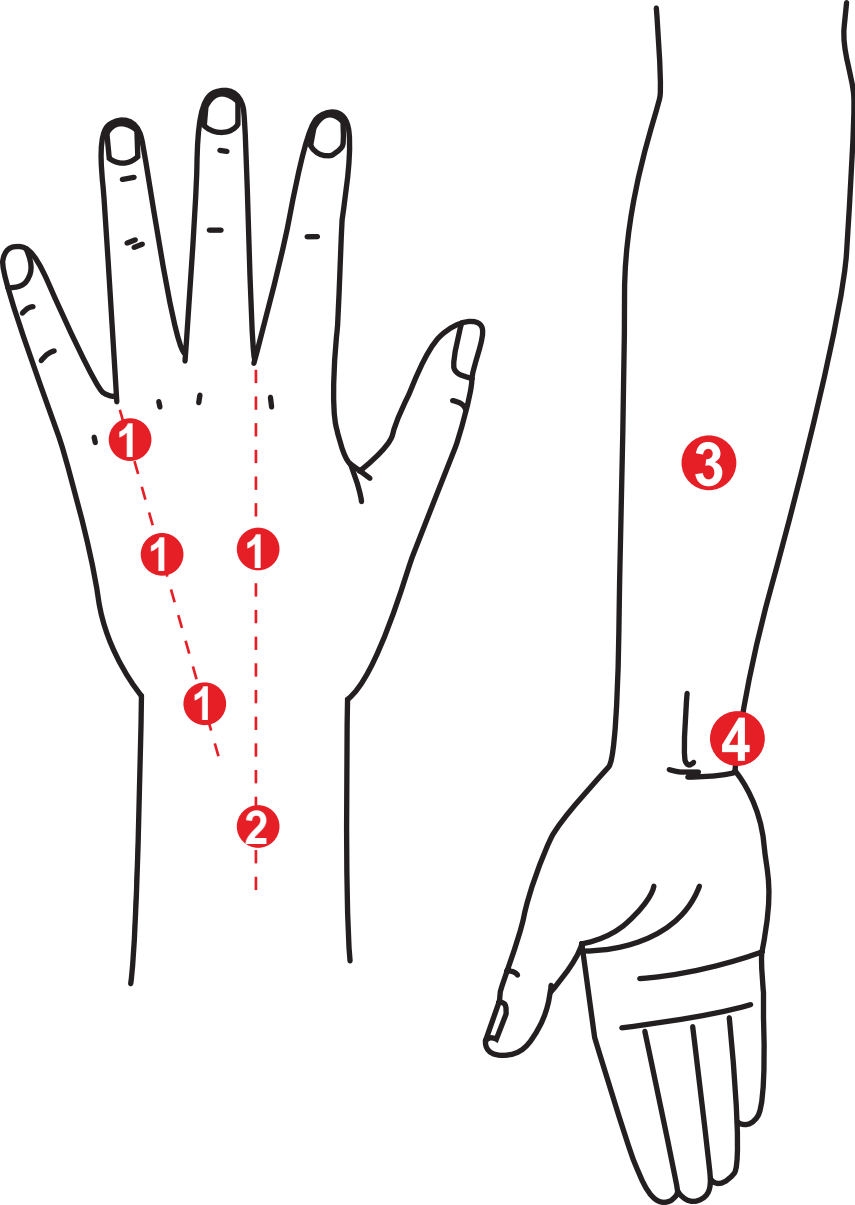

Place 2 way Jade Massage Heads
on the spots for **Fingers and Wrist.**

● spot for 2 way Jade Massage Heads


- 1: Pain of Fingers
- 2: Pain of Fingers, Fever
- 3: Improve Blood Circulation
- 4: Wrist Arthritis

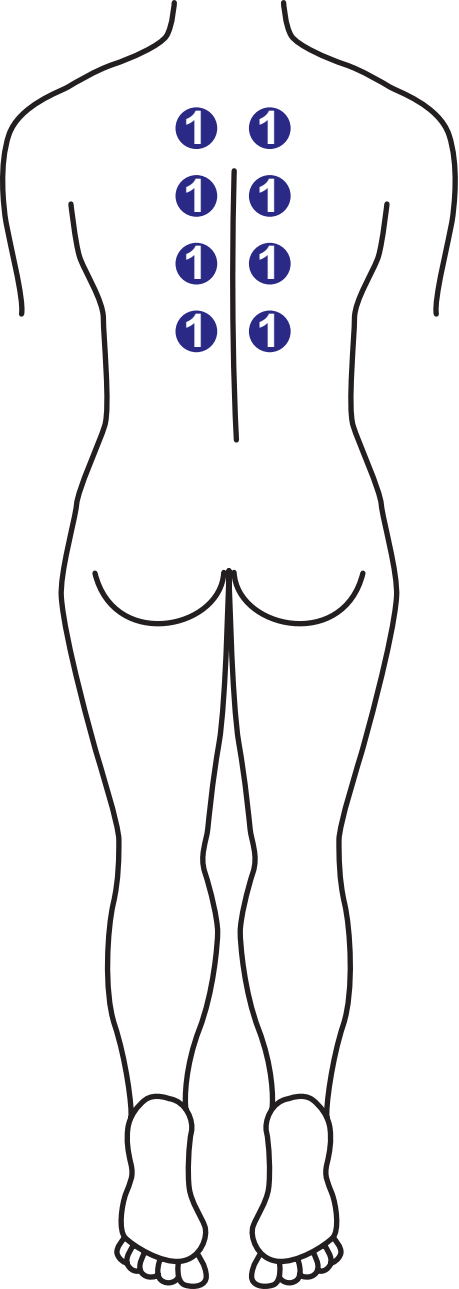



Place 5 way Jade Massage Heads
on the spots for **Aches and Pain.**

 spot for 5 way Jade Massage Heads

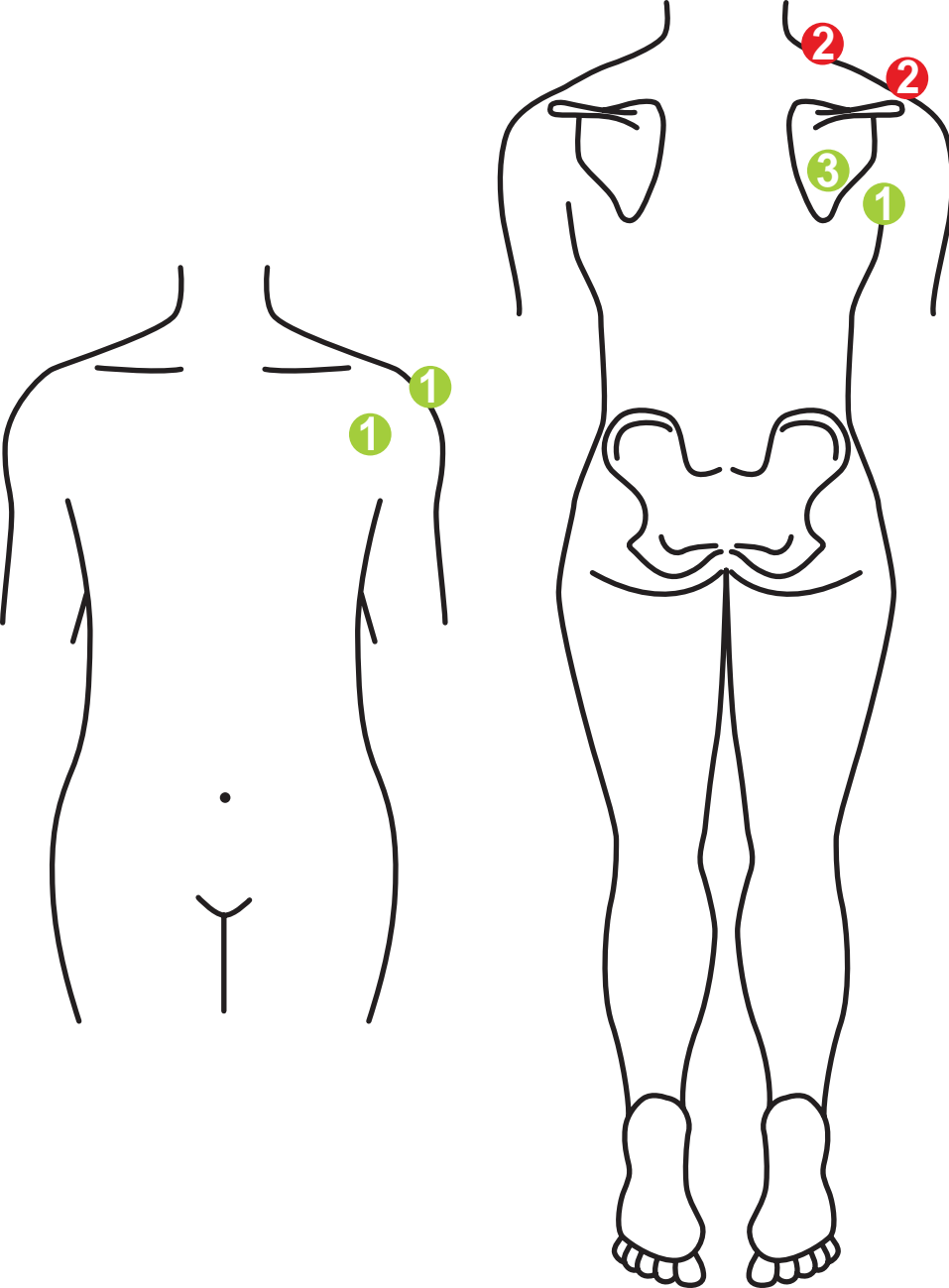
**1: Massage and Cauterize these
area several times with
5 way Jade Massage Heads.**



Place 2 & 7 way Jade Massage Heads
on the spots for **Shoulder Pain**.

-  spot for 2 way Jade Massage Heads
-  spot for 7 way Jade Massage Heads


- 1: For shoulder pain
- 2: For dizziness,
Neck muscle pain
- 3: For shoulder arthritis

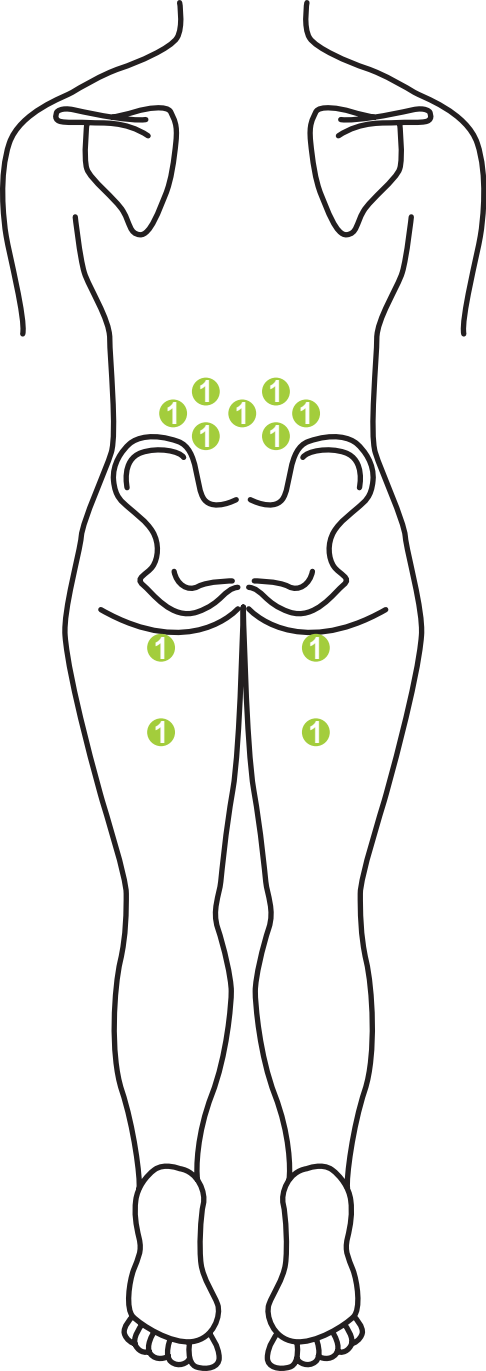

Place 7 way Jade Massage Heads
on the spots for **Lower Back Pain**.

 spot for 7 way Jade Massage Heads

1: Place 7 way Jade Massage Heads
on the floor and lie on top it.

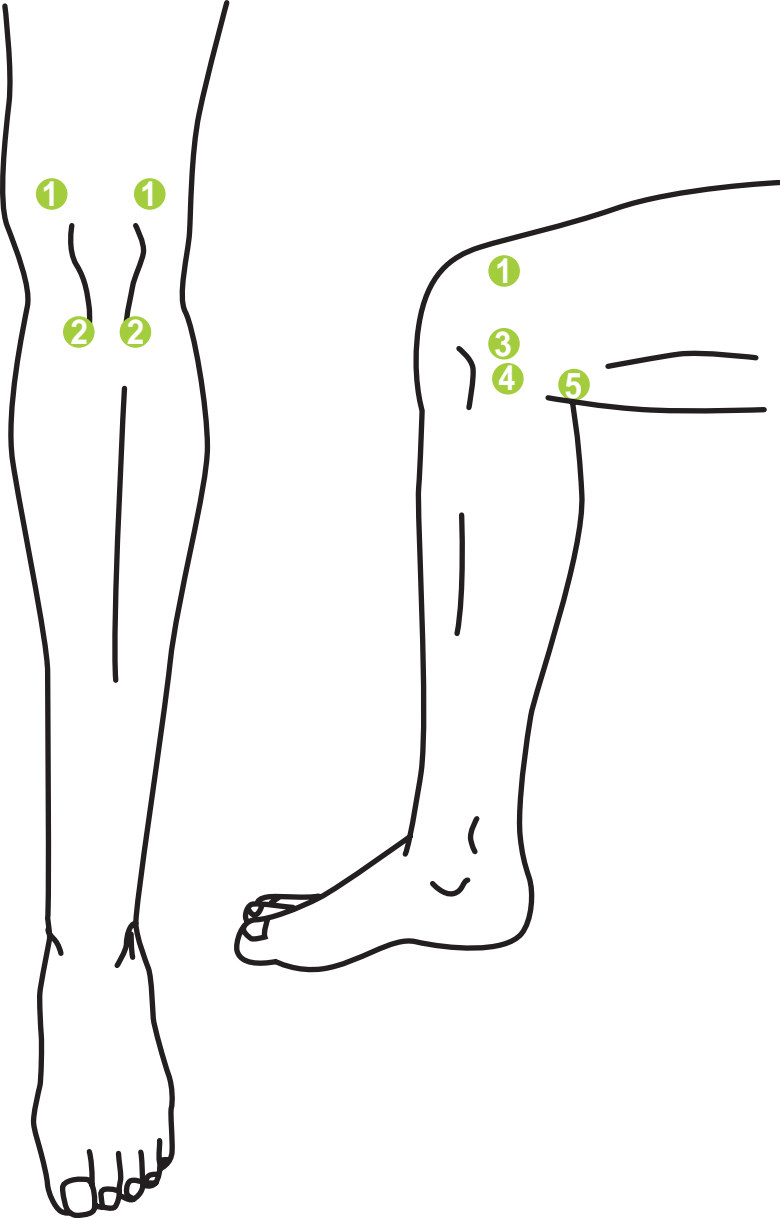

Place 7 way Jade Massage Heads
on the spots for **Knee Arthritis**.

 spot for 7 way Jade Massage Heads

- 1: Clear Blood
- 2: Knee Arthritis
- 3: Knee Problems
- 4: Knee Pain
- 5: Knee Arthritis, Leg Paralysis